



Smile
LINEO
will do it

LINEO IS A COMPLETE
ORTHODONTIC PACK:

ALIGNER +
SUPPORT IN THE DEVELOPMENT
OF THE TREATMENT PLAN,
ALL INCLUSIVE.

www.lineo.it

LINEO
THE HUMAN TOUCH ALIGNER
BY MICERIUMLAB

What's LINEO?

It is a **COMPLETE ORTODONTIC PACK:**

of the treatment plan

ALL INCLUSIVE

The aligner
alone isn't enough
to make your
patient happy.

Not even increasingly
sophisticated softwares
or intraoral scanners of
incredible precision, are.

You need support
of a clinical team,
experienced in invisible
orthodontics, which is able
to project and elaborate,
together with you, the right
treatment plan.

He is our clinical director: one of
the leading experts in Europe and
first in Italy for number of cases,
he coordinates the TEAM.
of our specialists in orthodontics

A practical guide on how to deal with the
most frequent problems of patients



LINEO transparent aligners are
ergonomic

Elegant aligner holder
included with patient's BOX

Why LINEO?

Because the patient's treatment plan will be elaborated according to your clinical needs by a TEAM OF SPECIALISTS coordinated by Dr. Aquilio, Clinical Manager and Scientific Director of LINEO.

The BOX that you will receive after accepting the treatment plan: inside it, you will find everything you need for the treatment of your patient.

The elegant BOX, with which, the first aligners of the treatment are delivered to the patient. The BOX contains the aligners, the guide to the correct use and the practical aligner holder.



A practical guide to the correct use of the aligner, included in the BOX to be delivered to patients



**MORE THAN
4000 CASES
HANDLED**

TRUST THE EXPERIENCE

Dr. Aquilio graduated in Dentistry and Dental Prosthetics at the University of Turin in 1992 and specialized in Orthodontics at the University of Cagliari in 2015. He attended for seven years the II Department of Orthodontics of the Hospital "Molinette" in Turin. He cooperates with the School of Specialization in Orthodontics of the University of Cagliari.

Is currently among the largest Europe's experts in invisible orthodontics with aligners



How does LINEO work?



IT IS SIMPLE AND THE
FIRST STEP **IS FREE**

Lineo is the ideal aligner for both orthodontists and non-orthodontic dentists. Orthodontists can count on a new aligner that is a valid alternative to other systems. Dentists can offer high quality orthodontic treatments to their patients.





“ I appreciate the dialogue between myself and the team of specialists with whom I developed my treatment plan: the feeling is of a clinical cooperation between professionals, not between a costumer and company. ”

Dr. Casu,
Dentist Specialist in Orthodontics, Sassari


Step 1 FREE TREATMENT , PLAN

The first Step is a FREE EVALUATION in which the CLINICAL TEAM LINEO provides information on the feasibility of your treatment plan, based on the photographs sent.

FREE STEP

Step 2 LINEOVISION

If you approve the PRELIMINARY PROJECT on the case, the CLINICAL TEAM develops LINEO VISION, the virtual care plan based on the patient's impressions.

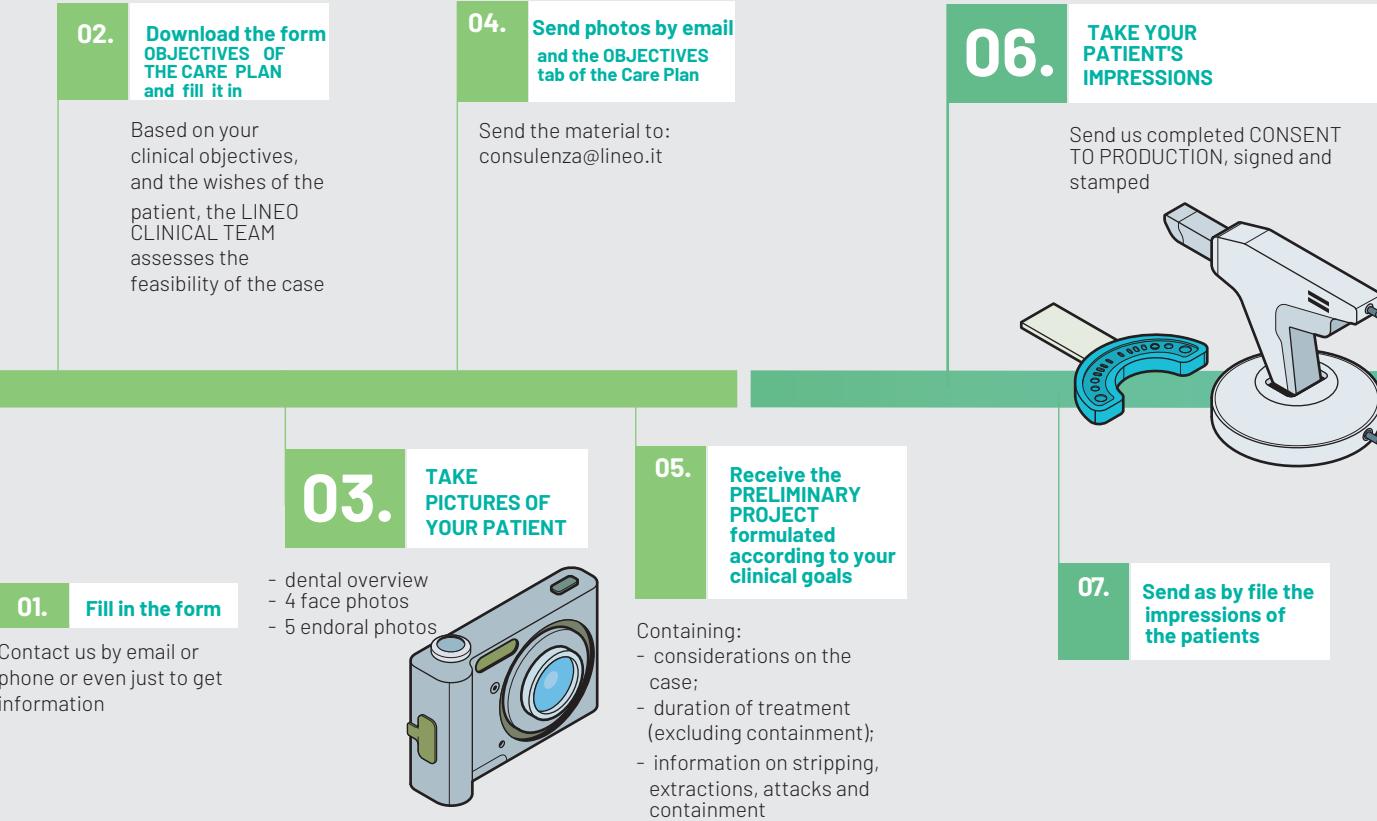




Step 3 TREATMENT

If you accept your LINEOVISION, we proceed with the production of the aligners and the shipment, directly to your Practice.

PAID PHASES

08. CASE DEVELOPMENT

The LINEO CLINICAL TEAM develops LINEOVISION, it is the virtual care plan for the your case

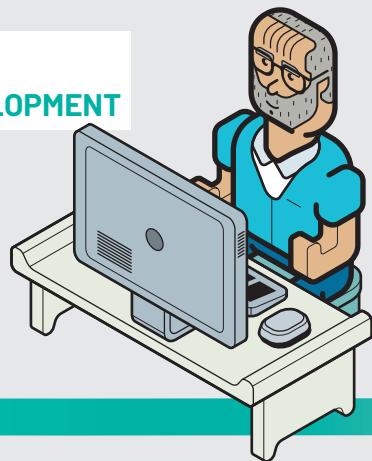

10. Accept LINEOVISION

If you accept your LINEOVISION, the production of the aligners proceeds and the following shipment is made directly to your practice

12. Contact us when you want

We are available as many times as you wish, call us or send us a photo by email

09.

RECEIVE YOUR LINEOVISION, YOUR VIRTUAL CARE PLAN

It is a three-dimensional and interactive animation to be approved and shown to the patient



11.

Send photos every 8 aligners

Every 8 aligners it is recommended to send 5 intraoral photos according to the usual projections without aligners and 2 intraoral photos with aligners

The advantages FOR THE DOCTOR

A QUALITY CARE PLAN

Support in the development of the care plan is included, thanks to a coordinated CLINICAL TEAM by one of the leading experts in Europe and first in Italy for the number of cases handled.



DIRECT ASSISTANCE

We will be available as often as you wish to call us; be it for a clinical comparison of the treatment, a clarification or a simple question.

THE EXCLUSIVE CONTACT WITH YOUR PATIENT

Lineo will never come into contact with your patient: you will be the sole and only referent of the treatment.



The advantages FOR THE PATIENT

AN ERGONOMIC SHAPE OF THE ALIGNER

The LINEO aligner is characterized, in the vestibular part, by a gingival festoon for impeccable aesthetics.

On the contrary, the edge of the lingual/palatine part, slides about 3 mm from the gingival festoon, increasing the strength of the aligner and the patients comfort.




QUALITY TREATMENT

You can offer your patients a treatment developed personally by one of the leading experts in Europe and the first in Italy for the number of cases treated.

THE DENTIST YOU TRUST

The patient can benefit from the orthodontic treatment with aligners at his or hers own dentist of confidence.

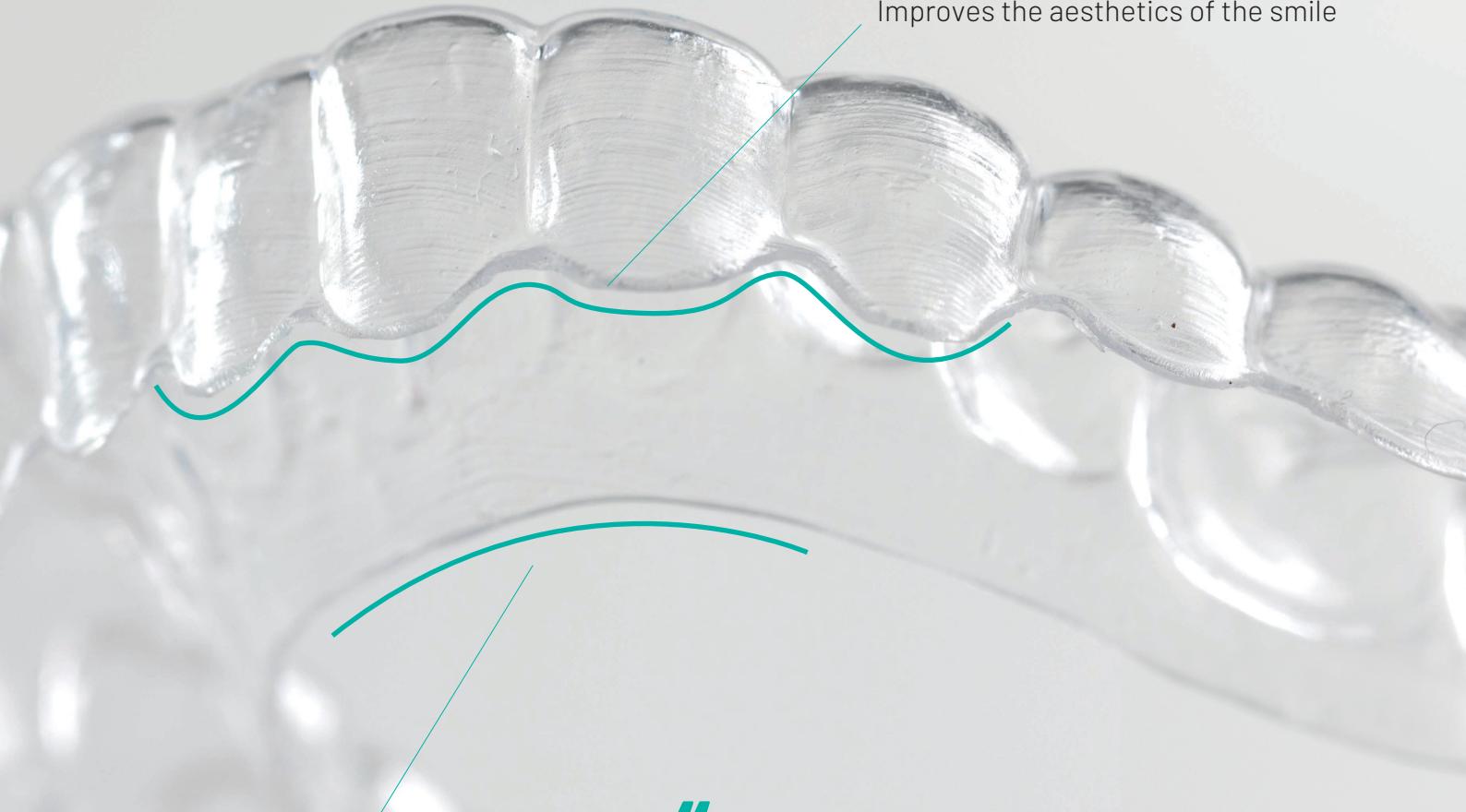

CURVED LANGUAGE EDGE

Increases robustness and comfort of the patient

VESTIBULAR EDGE

Improves the aesthetics of the smile

“

I am in the middle of the treatment and I am very happy with how my smile is changing: the aligners are transparent, comfortable and very easy to wear.

”

Giovanna D. 23 years old, Turin

Auxiliary Elements

LINEO

TO SOLVE EVERY KIND OF CASE

SIMPLE AND COMPLEX CASES

HORIZONTAL CHAMFERED ATTACHMENT

INDICATIONS: facilitates the extrusion movement.
Specifications for upper and lower incisors.

BEVELLED VERTICAL PYRAMIDAL ATTACHMENT


INDICATIONS: rotations. It is used in pairs: one on the surface of the vest and one on the surface of the blade / lingual to create a pair of forces. Specifications for upper and lower molars.

VERTICAL RECTANGULAR ROUNDED ATTACHMENT

INDICATIONS: complex/changing rotations / increasing aligner retention. Used on molars, premolars, upper and lower canines.

CUTS FOR ELASTIC TRACTION

INDICATIONS: used to carry out elastic traction of 2nd and 3rd class.
These are incisions made on the aligners at the level of canines and upper and lower molars.

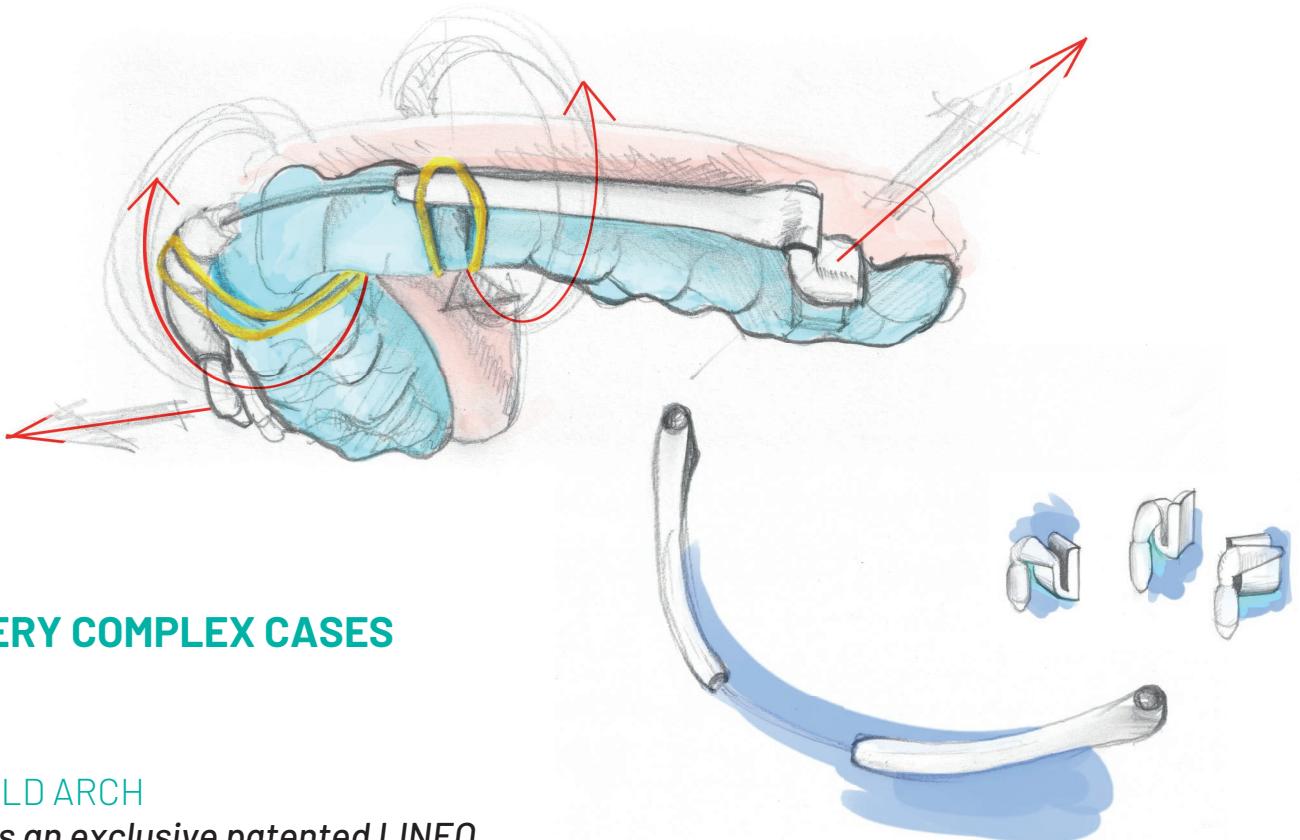

VERY COMPLEX CASES


WILD ARCH

It is an exclusive patented LINEO

INDICATIONS: posterior cross-bites, severely contracted arches and severe deep bites.

Customized removable device:
is used at night.
Easy to handle for the patient.

Thanks to its beta-titanium core,
it can increase by 100%, the
performance of the aligners, with
regard to arch expansions and the
intrusion of upper and lower canines/
incisors.

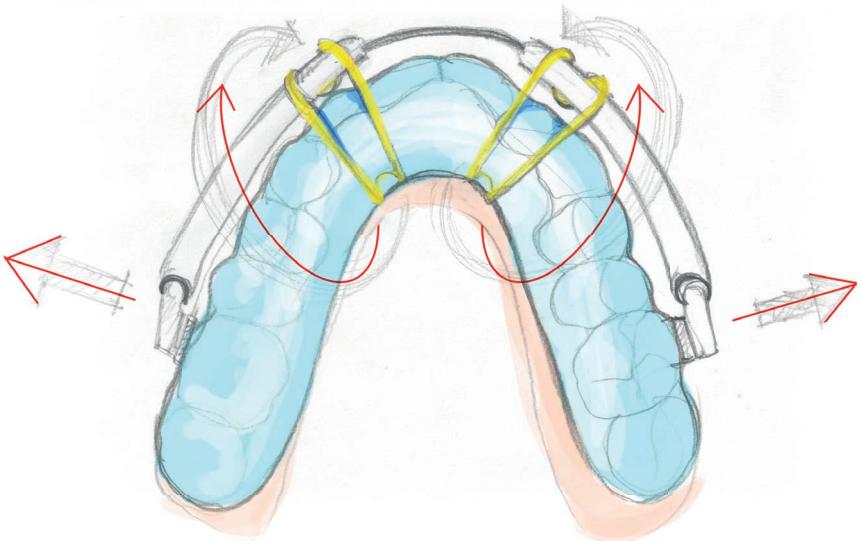

Become a LINEO PARTNER

Give your patients the opportunity
to improve their smile.

With minimal impact in terms of time and know-how on your part, you can turn LINEO into a marketing tool to increase the number of your patients, offering them a new high quality aligner at a very affordable price.

Call us at the number

Staff Clinico LINEO: T +39 011 044.67.47

Business Care LINEO: T +39 0185 79.121

Or go to

www.LINEO.IT

www.miceriumlab.com

“Compared to other aligners, I have more freedom on the number of annual cases: whether they are many or few, the price is always very affordable and very high quality.”

Dr. Terrando, Dentist, Turin

The elegant BOX with which the first aligners of the treatment are delivered to the patient. The BOX contains the aligners, the guide to the correct use and the practical aligner holder.



LINEO

THE HUMAN TOUCH ALIGNER
BY MICERIUM LAB



www.lineo.it